



# EDGEWOOD CHEMICAL BIOLOGICAL CENTER

U.S. ARMY RESEARCH, DEVELOPMENT AND ENGINEERING COMMAND  
Aberdeen Proving Ground, MD 21010-5424

ECBC-TR-1450

## QUANTIFICATION OF TETRAMETHYLENEDISULFOTETRAMINE (TETS) IN VARIOUS FOOD MATRICES BY SOLID-PHASE EXTRACTION LIQUID CHROMATOGRAPHY-ION TRAP MASS SPECTROMETRY

Sue Y. Bae  
Li Kong

RESEARCH AND TECHNOLOGY DIRECTORATE

April 2017

Approved for public release: distribution unlimited.



#### Disclaimer

The findings in this report are not to be construed as an official Department of the Army position unless so designated by other authorizing documents.

REPORT DOCUMENTATION PAGE				Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 h per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Department of Defense, Washington Headquarters Services, Directorate for Information Operations and Reports (0704-0188), 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to any penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number. PLEASE DO NOT RETURN YOUR FORM TO THE ABOVE ADDRESS.					
1. REPORT DATE (DD-MM-YYYY) XX-04-2017		2. REPORT TYPE Final		3. DATES COVERED (From - To) Mar 2012 – Jun 2012	
4. TITLE AND SUBTITLE Quantification of Tetramethylenedisulfotetramine (TETS) in Various Food Matrices by Solid-Phase Extraction Liquid Chromatography–Ion Trap Mass Spectrometry				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S) Bae, Sue Y. and Kong, Li				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Director, ECBC, ATTN: RDCB-DRC-C, APG, MD 21010-5424				8. PERFORMING ORGANIZATION REPORT NUMBER ECBC-TR-1450	
9. SPONSORING / MONITORING AGENCY NAME(S) AND ADDRESS(ES) Defense Threat Reduction Agency, 8725 John J. Kingman Road, MSC 6201, Fort Belvoir, VA 22060-6201				10. SPONSOR/MONITOR'S ACRONYM(S) DTRA	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION / AVAILABILITY STATEMENT Approved for public release: distribution unlimited.					
13. SUPPLEMENTARY NOTES					
14. ABSTRACT: Tetramethylenedisulfotetramine (TETS), commonly known as tetramine, is a highly neurotoxic rodenticide (the LD <sub>50</sub> , or lethal dose for 50% of test subjects, is 0.1 mg/kg). TETS has been used in hundreds of deliberate and accidental food poisoning events in China. In this report, we describe an extraction method for the quantitation of TETS as spiked into various food matrices, including fruit juices, egg, hot dog, chicken nuggets, turkey deli meat, and ground beef. Quantitation by liquid chromatography–ion trap mass spectrometry was based upon selected reaction monitoring of mass-to-charge ratio ( $m/z$ ) 347 → 268 and consecutive reaction monitoring of $m/z$ 347→268→175. An ~5 mg TETS concentration was spiked into each food sample. The linear range of quantitation for TETS was 0.5–20 µg/mL. Total recoveries (and percent relative standard deviations) for TETS in various food samples are reported.					
15. SUBJECT TERMS Tetramethylenedisulfotetramine (TETS)      Food Selected reaction monitoring (SRM)      Liquid chromatography–ion trap mass spectrometry (LC–IT-MS) Consecutive reaction monitoring (CRM)					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT  UU	18. NUMBER OF PAGES  38	19a. NAME OF RESPONSIBLE PERSON Renu B. Rastogi
a. REPORT U	b. ABSTRACT U	c. THIS PAGE U			19b. TELEPHONE NUMBER (include area code) (410) 436-7545

Blank

## **PREFACE**

This work was started in March 2012 and completed in June 2012.

The use of either trade or manufacturers' names in this report does not constitute an official endorsement of any commercial products. This report may not be cited for purposes of advertisement.

This report has been approved for public release.

## **Acknowledgments**

We thank Kim Green (Director, Office of the Chief Scientist, U.S. Department of Agriculture [USDA]; Washington, DC) and Nick Bauer (USDA Food Safety and Inspection Service) for their initial funding of this project.

Blank

## CONTENTS

1.	INTRODUCTION .....	1
2.	EXPERIMENTAL METHODS.....	2
2.1.	Reagents and Chemicals .....	2
2.2.	Instrumentation .....	2
2.3.	Sample Preparation and Extraction Procedure .....	3
2.3.1	Juice Sample Preparation.....	3
2.3.2	Egg Sample Preparation.....	3
2.3.3	Hot Dog, Chicken Nuggets, Turkey Deli Meat, and 80/20 Ground Beef Preparation .....	4
2.3.4	Extraction Procedures .....	4
3.	RESULTS AND DISCUSSION .....	5
3.1.	Apple and Orange Juices .....	5
3.2.	Egg Beaters Egg Whites; Whole Egg; and USDA Egg White, Salted Egg Yolk 2551, and P. Sugared Egg Yolk C24410.....	7
3.3.	Hot Dog, Chicken Nuggets, Turkey Deli Meat, and 80/20 Ground Beef.....	9
4.	CONCLUSION.....	11
	LITERATURE CITED .....	13
	ACRONYMS AND ABBREVIATIONS .....	15
	APPENDIX: SELECTED REACTION MONITORING (SRM) AND CONSECUTIVE REACTION MONITORING (CRM) CHROMATOGRAMS AND MASS SPECTRA OF TETRAMETHYLENEDISULFOTETRAMINE (TETS) EXTRACTED FROM VARIOUS FOOD MATRICES .....	17

## FIGURES

1.	Chemical structure of TETS .....	2
2.	Schematic of RediSep Rf normal-phase silica column.....	4
3.	External calibration curve for TETS in acetonitrile and water for juice analyses .....	5
4.	Representative (a) SRM and (b) MS results for TETS extracted from apple juice matrix.....	6
5.	External calibration curve for TETS in acetonitrile and water for various egg sample analyses .....	7
6.	Representative (a) CRM and (b) MS results for TETS extracted from Egg Beaters egg whites matrix.....	8
7.	External calibration curve for TETS in acetonitrile and water for various hot dog, chicken nugget, turkey deli meat, and 80/20 ground beef sample analyses .....	9
8.	Representative (a) CRM and (b) MS for TETS extracted from 80/20 ground beef matrix.....	10

## TABLES

1.	LC–IT-MS Parameters .....	3
2.	Percent Recovery of TETS Extracted from Juice Matrices .....	7
3.	Percent Recovery of TETS Extracted from Various Egg Matrices .....	9
4.	Percent Recovery of Extracted TETS from Hot Dog, Chicken Nugget, Turkey Deli Meat, and 80/20 Ground Beef Matrices .....	11

# QUANTIFICATION OF TETRAMETHYLENEDISULFOTETRAMINE (TETS) IN VARIOUS FOOD MATRICES BY SOLID-PHASE EXTRACTION LIQUID CHROMATOGRAPHY–ION TRAP MASS SPECTROMETRY

## 1. INTRODUCTION

Tetramethylenedisulfotetramine (TETS), commonly known as tetramine, is a highly neurotoxic rodenticide (the LD<sub>50</sub>, or lethal dose for 50% of test subjects, is 0.1 mg/kg).<sup>1</sup> TETS has been used in hundreds of deliberate and accidental food poisoning events in China. Banned rodenticides illicitly imported into the United States have been responsible for numerous intentional and unintentional human poisonings.<sup>2–4</sup> TETS is slightly soluble in water (0.25 mg/mL),<sup>5</sup> and it is dangerous when used to deliberately contaminate food or water. TETS is not absorbed through the skin; the most common route of exposure is ingestion of contaminated foods. Thus, the development of a reliable extraction and detection technique for TETS in different foods is essential: when accidental and intentional poisonings occur, responders must be able to identify the presence of potentially toxic substances in foods.

Most of the TETS detection in food and biological sample matrices has been achieved using gas chromatography coupled with mass spectrometry (GC–MS),<sup>6–10</sup> and limited literature exists regarding the direct detection of TETS in food via liquid chromatography coupled with mass spectrometry (LC–MS).<sup>11</sup> Several reports have been published describing sample preparation and cleanup techniques that rely on solid-phase microextraction,<sup>9,12</sup> however, the longer extraction time (>60 min) and analysis time (>10 min) make this technique inefficient for processing large numbers of samples.

This report documents the efforts of the Agent Chemistry Team from the Research and Technology Directorate of the U.S. Army Edgewood Chemical Biological Center (ECBC; Aberdeen Proving Ground, MD) in developing new extraction and analytical detection methodologies using liquid chromatography–ion trap mass spectrometry (LC–IT-MS). The objective of this task was to provide development and laboratory support for extraction of TETS (Figure 1) from various food samples. This included detection and quantitative and qualitative analyses of complex matrices, such as foods with high salt and fat contents. In support of this objective, we examined 11 food samples using individual agents, including apple and orange juices; Egg Beaters processed egg whites; U.S. Department of Agriculture (USDA) egg white, salted egg yolk 2551, and p. sugared egg yolk C24410; whole egg; chicken nuggets; hot dog; precooked, 99% fat free turkey deli meat; and 80/20 hamburger meat (80% lean and 20% fat). The choice of food types arose from collaboration and conversations with USDA personnel and represent items that are commonly associated with school lunch programs. Foods were cleaned using commercially available normal-phase separation columns.

The use of LC–IT-MS, or any comparable high-resolution tandem mass spectrometry, has become common. From an affordability standpoint, this technology is currently within reach for most laboratories. For this work, extracted TETS was analyzed using LC–IT-MS, and percent recoveries were calculated from external calibration curves.

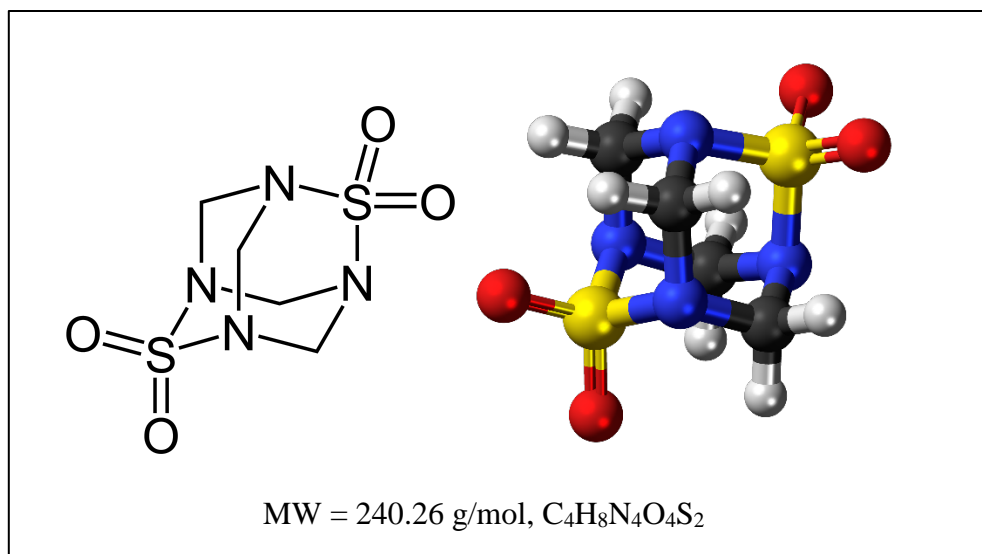


Figure 1. Chemical structure of TETS.

## 2. EXPERIMENTAL METHODS

### 2.1 Reagents and Chemicals

The TETS (>98% purity) was provided by ECBC. All reagents and solvents were LC–MS grade. Acetonitrile, water, and triethylamine (TEA) were purchased from Sigma-Aldrich (St. Louis, MO). A RediSep Rf normal-phase silica gel column (5 g) was obtained from Teledyne Isco (Lincoln, NE). Juicy Juice apple juice (Harvest Hill Beverage Company; Stamford, CT), Minute Maid orange juice (Coca Cola Company; Atlanta, GA), Egg Beaters processed egg whites (ConAgra Foods; Chicago, IL), whole egg, Smart Option chicken nuggets (Food Lion; Salisbury, NC), Esskay Oriole hot dog (Smithfield Foods; Smithfield, VA), Buddig original thin turkey deli meat (Carl Buddig and Company; Homewood, IL), and 80/20 ground beef food samples were purchased from a local grocery store (Food Lion; Edgewood, MD). Four egg samples were provided from USDA: egg white, salted egg yolk 2551, and p. sugared egg yolk C24410.

### 2.2 Instrumentation

All samples were characterized using an LTQ ion trap mass spectrometer (Thermo Fisher Scientific; Waltham, MA) equipped with an electrospray ionization (ESI) interface. The operating parameters are tabulated in Table 1.

Table 1. LC–IT-MS Parameters

Column: Phenomenex Gemini-NX 3 μm, C18 50 × 2 mm				
Mobile phase A: water				
Mobile phase B: acetonitrile				
Gradient	Time (min)	A (%)	B (%)	Flow Rate (μL/min)
	0	75	25	350
	1	75	25	
	3	5	95	
	4.5	5	95	
	4.6	75	25	
	7	75	25	
Injection volume: 5 μL				
Sheath gas flow rate: 50 L/min				
Auxiliary gas flow rate: 20 L/min				
Sweep gas flow rate: 10 L/min				
Ion spray voltage: 3.5 kV				
Capillary temperature: 300 °C				
Capillary voltage: −50 V				
Tube lens voltage: −86 V				
Ionization mode: ESI				

## 2.3 Sample Preparation and Extraction Procedure

### 2.3.1 Juice Sample Preparation

Juicy Juice apple juice and Minute Maid orange juice were purchased from a local grocery store. Two milliliters of each sample was spiked with ~5 mg of TETS and diluted with 20 mL of acetonitrile. The food sample was then passed through a RediSep Rf column, and the eluents were collected.

### 2.3.2 Egg Sample Preparation

Original Egg Beaters egg whites and whole eggs were purchased from a local grocery store. USDA egg white, salted egg yolk 2551, and p. sugared egg yolk C24410 were provided by USDA. Approximately 5 g of each food sample was spiked with the desired TETS quantity and diluted with 10 mL of acetonitrile. The mixture was centrifuged for 15 min at 6000 rpm, and the supernatant was decanted. A second 10 mL portion of acetonitrile was added, and the mixture was vortexed or sonicated for 1 min and again centrifuged for 15 min at 6000 rpm. The supernatant was removed, combined, and passed through a RediSep Rf column. The eluents were collected for analysis.

### 2.3.3 Hot Dog, Chicken Nuggets, Turkey Deli Meat, and 80/20 Ground Beef Preparation

Esskay hot dogs, Smart Option chicken nuggets, Buddig turkey deli meat, and 80/20 ground beef were purchased from a local grocery store. Approximately 5 g of each food sample was spiked with TETS and diluted with 10 mL of acetonitrile. The whole sample was homogenized using a Kinematica Polytron homogenizer (Kinematica AG; Luzern, Switzerland) at 20,000 rpm for 1–2 min. The mixture was then centrifuged for 15 min at 6000 rpm, and the supernatant was removed. A second 10 mL portion of acetonitrile was added, and the sample was vortexed or sonicated for 1 min and centrifuged for 15 min at 6000 rpm, and the supernatant was removed. The combined supernatant was passed through a RediSep Rf column, and the eluents were collected for analysis.

### 2.3.4 Extraction Procedures

A packed RediSep Rf normal-phase silica gel column (shown in Figure 2) was used in this study to separate TETS from the food samples. First, the column was eluted with 25 mL of 0.1% TEA/acetonitrile, using in-house air to pass the solution through the column, and the 0.1% TEA/acetonitrile solution was collected for later use. Second, the supernatant was passed through the column, and the sample was collected. Third, 1 mL of 0.1% TEA/acetonitrile solution was added to the column and pushed slightly into the silica gel until 1 mL of the solution just cleared the top of the silica gel. This step was repeated four times. Finally, the remaining 0.1% TEA/acetonitrile solution was added to the column and passed through the bed. A small aliquot of the extracted solution was filtered through a 0.45  $\mu\text{m}$  poly(tetrafluoroethylene) membrane filter. The final solution was diluted with mobile phase and analyzed using LC–MS.



Figure 2. Schematic of RediSep Rf normal-phase silica column.

### 3. RESULTS AND DISCUSSION

#### 3.1 Apple and Orange Juices

The quantitation was evaluated by linear regression with an external calibration curve (Figure 3) that ranged from 0.10 to 20  $\mu\text{g/mL}$  with  $n \geq 9$  measurements per standard. A standard curve was prepared at the beginning of each sequence run. In the presence of heat (i.e., mass spectrometer collision cell), TETS forms the TETS dimer,<sup>11</sup> which in turn becomes protonated before losing methylene to form an ion of mass-to-charge ratio ( $m/z$ ) 347. The possible product ions were  $m/z$  268, 227, 175, 148, and 134. TETS was analyzed using selected reaction monitoring (SRM) of the  $m/z$  347  $\rightarrow$  268 transition. A representative SRM for TETS in an apple juice sample is shown in Figure 4a. The mass spectrum, shown in Figure 4b, exhibits mass ions at  $m/z$  268.09 due to  $[\text{M}-\text{SO}_2\text{NH}]^-$  for TETS. The percent recovery was calculated based on an external calibration curve for TETS, and the results revealed a >95% recovery of TETS from apple juice and orange juice (Table 2).

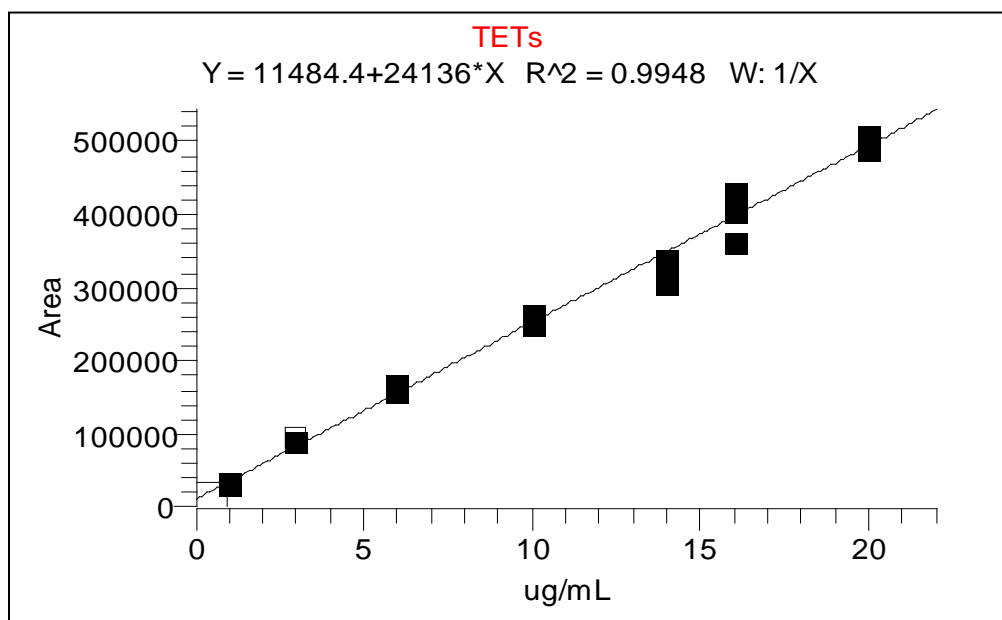


Figure 3. External calibration curve for TETS in acetonitrile and water for juice analyses.

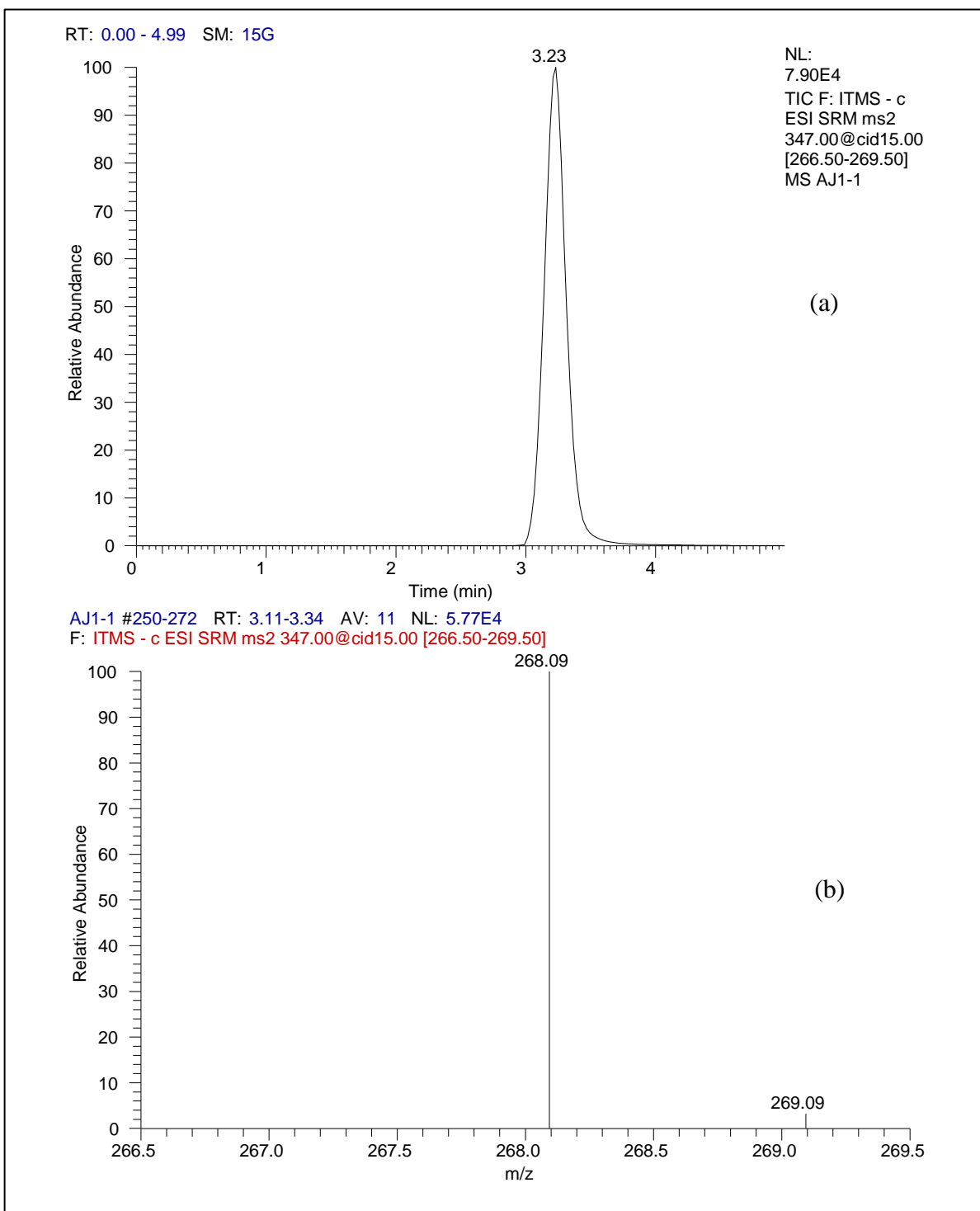


Figure 4. Representative (a) SRM and (b) MS results for TETS extracted from apple juice matrix.

Table 2. Percent Recovery of TETS Extracted from Juice Matrices

Sample	pH	Recovery (%)
Apple juice	3.5	100.7 ± 2.3
Orange juice	4.0	96.9 ± 2.0

### 3.2 Egg Beaters Egg Whites; Whole Egg; and USDA Egg White, Salted Egg Yolk 2551, and P. Sugared Egg Yolk C24410

The quantitation was evaluated by linear regression with an external calibration curve (Figure 5) that ranged from 0.10 to 20 µg/mL with  $n \geq 9$  measurements per standard. A standard curve was prepared at the beginning of each sequence run. TETS was analyzed using consecutive reaction monitoring (CRM), and the  $m/z$  347 → 268 → 175 transitions were monitored. A representative CRM for TETS in an Egg Beaters sample is shown in Figure 6a. The mass spectrum, shown in Figure 6b, exhibits mass ions at  $m/z$  175.05 due to  $[M-CH_4N_2O_4S_2]^-$  for TETS. The percent recovery was calculated based on an external calibration curve for TETS, and the results showed >88 % recovery of TETS from Egg Beaters egg whites, whole egg, and USDA egg white. However, the percent recoveries of TETS from USDA salted egg yolk 2551 and p. sugared egg yolk C24410 were <90% (Table 3). Two of the USDA egg yolk samples were very thick and sticky. A procedure that includes multiple centrifuge and extraction steps could improve the percent recovery of TETS from the salted and p. sugared egg yolk samples.

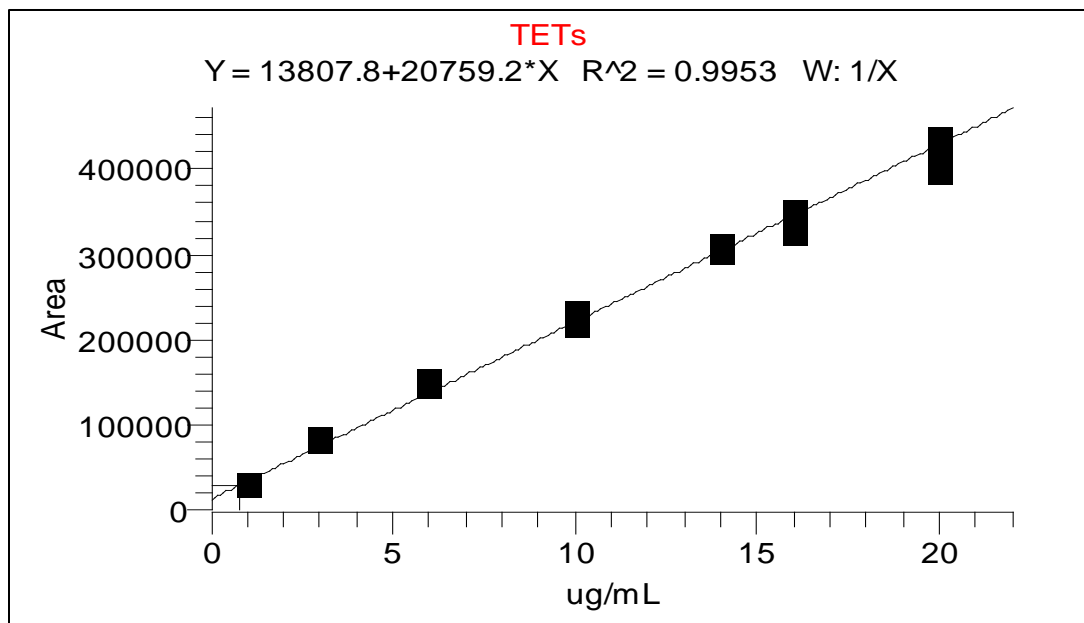


Figure 5. External calibration curve for TETS in acetonitrile and water for various egg sample analyses.

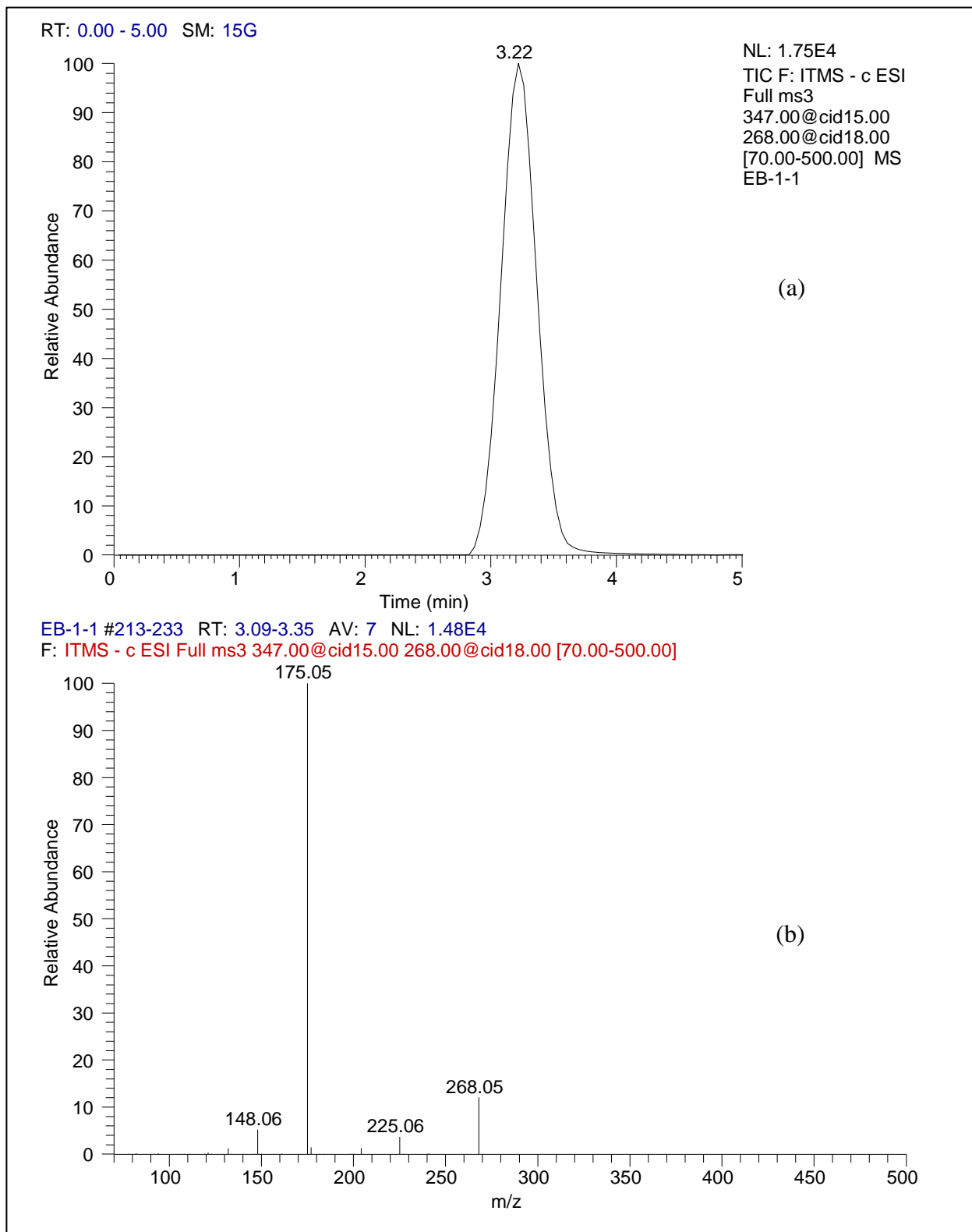


Figure 6. Representative (a) CRM and (b) MS results for TETS extracted from Egg Beaters egg whites matrix.

Table 3. Percent Recovery of TETS Extracted from Various Egg Matrices

Sample	pH	Recovery (%)
Whole egg	7	98.7 ± 1.4
USDA egg white	7	96.1 ± 2.7
Egg Beaters egg white	7	91.5 ± 1.6
USDA p. sugared egg yolk C24410	7	88.9 ± 3.2
USDA salted egg yolk 2551	7	79.1 ± 8.4

### 3.3 Hot Dog, Chicken Nuggets, Turkey Deli Meat, and 80/20 Ground Beef

The quantitation was evaluated by linear regression with an external calibration curve (Figure 7) that ranged from 0.10 to 20 µg/mL with  $n \geq 9$  measurements per standard. A standard curve was prepared at the beginning of each sequence run. TETS was analyzed using CRM, and the  $m/z$  347 → 268 → 175 transitions were monitored. A representative CRM for TETS in an 80/20 ground beef sample is shown in Figure 8a. The mass spectrum, shown in Figure 8b, exhibits mass ions at  $m/z$  175.05 due to  $[M-CH_4N_2O_4S_2]^-$  for TETS. Percent recovery was calculated based on an external calibration curve for TETS, and results showed that >88% recovery of TETS was realized from hot dog, chicken nuggets, turkey deli meat, and 80/20 ground beef (Table 4).

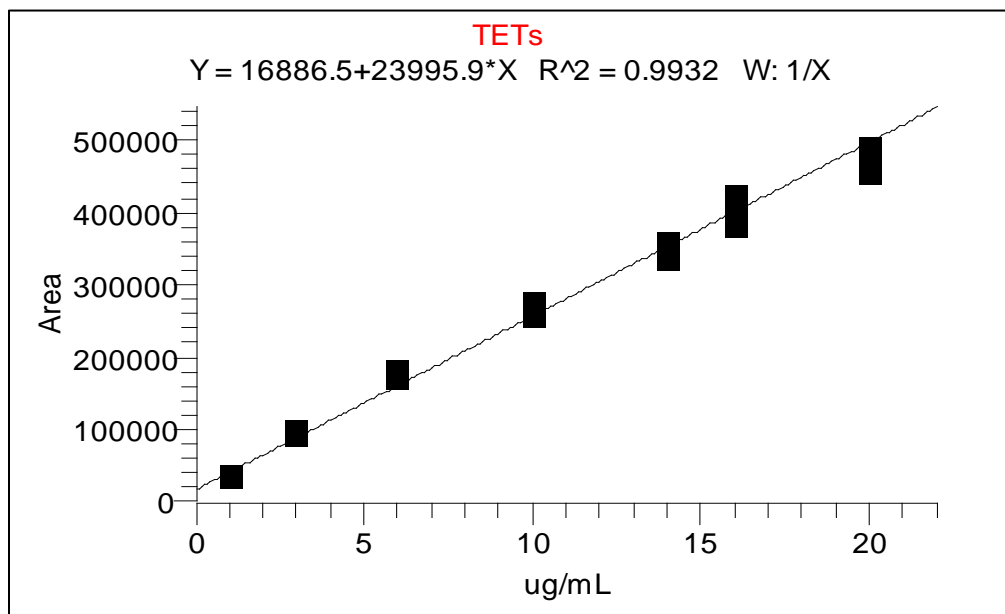


Figure 7. External calibration curve for TETS in acetonitrile and water for various hot dog, chicken nugget, turkey deli meat, and 80/20 ground beef sample analyses.

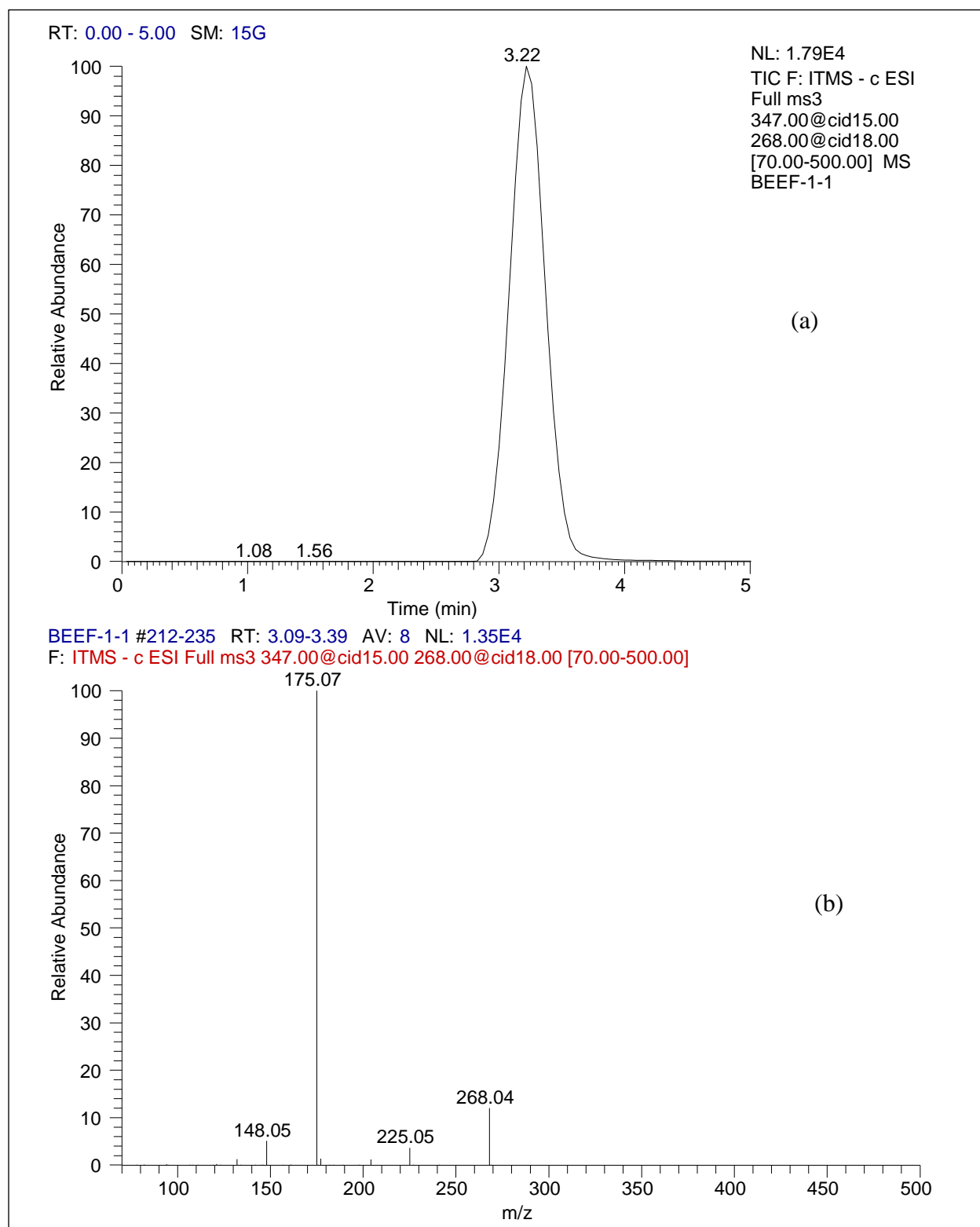


Figure 8. Representative (a) CRM and (b) MS for TETS extracted from 80/20 ground beef matrix.

Table 4. Percent Recovery of Extracted TETS from Hot Dog, Chicken Nugget, Turkey Deli Meat, and 80/20 Ground Beef Matrices

Sample	pH	Recovery (%)
Hot dog	7	93.5 $\pm$ 1.6
Chicken nuggets	7	87.4 $\pm$ 1.3
Turkey deli meat	7	90.1 $\pm$ 1.2
80/20 ground beef	7	88.8 $\pm$ 3.3

#### 4. CONCLUSION

The techniques developed for extracting TETS from various food matrices were accomplished, and recoveries >95% were obtained from most of the liquid and semi-liquid food matrices. The recoveries were >88% from all of the solid food matrices, including both of the USDA egg yolk samples (C24410 and 2551). We should be able to optimize the cleanup and extraction procedures to achieve greater TETS recovery from the high-fat, -salt, and -sugar food matrices. At present, we are working on the optimization process.

The extraction method was easy to use, and from it, we could determine the amounts of TETS in various complex food matrices. TETS is a relatively persistent environmental contaminant due to its high stability in water. This extraction method is applicable to a wide variety of complex foods, yet it requires minimal sample handling. The method is simple and inexpensive; thus, food samples could be extracted in the field by individuals who have minimal training. This would be ideal for quickly identifying potentially toxic food products anywhere in the United States.

Blank

## LITERATURE CITED

1. Croddy, E. Rat Poison and Food Security in the People's Republic of China: Focus on Tetramethylene Disulfotetramine (Tetramine). *Arch. Toxicol.* **2004**, 78 (1), 1–6.
2. Poisoning by an Illegally Imported Chinese Rodenticide Containing Tetramethylenedisulfotetramine—New York City, 2002. *MMWR Morb. Mortal Wkly Rep.* **2003**, 52 (10), 199–201.
3. Poisonings Associated with Illegal Use of Aldicarb as a Rodenticide—New York City, 1994–1997. *MMWR Morb. Mortal Wkly Rep.* **1997**, 46 (41), 961–963.
4. Nelson, L.S.; Perrone, J.; DeRoos, F.; Stork, C.; Hoffman, R.S. Aldicarb Poisoning by an Illicit Rodenticide Imported into the United States: Tres Pasitos. *J. Toxicol. Clin. Toxicol.* **2001**, 39 (5), 447–452.
5. Banks, C.N.; Yang, D.; Lein, P.J.; Rogawski, M.A. Tetramethylenedisulfotetramine. In *Encyclopedia of Toxicology*, 3rd ed.; Elsevier: New York, 2014; Vol. 4, pp 509–511.
6. Rong, W.G.; Liu, H.-I.; Chen, B.; Zhu, F.; Ji, W.-I.; Ma, Y.J. Determination of Four Rodenticides in Foods by Gas Chromatography–Tandem Mass Spectrometry. *Fenxi Shiyanshi [Chinese Journal of Analysis Laboratory]* **2013**, 32 (1), 73–77.
7. De Jager, L.S.; Perfetti, G.A.; Diachenko, G.W. *J. Sep. Sci.* **2009**, 32, 1081.
8. De Jager, L.S.; Perfetti, G.A.; Diachenko, G.W. Stir Bar Sorptive Extraction-Gas Chromatography-Mass Spectrometry Analysis of Tetramethylene Disulfotetramine in Food: Method Development and Comparison to Solid-Phase Microextraction. *Anal. Chim. Acta* **2009**, 635, 162–166.
9. De Jager, L.S.; Perfetti, G.A.; Diachenko, G.W. Analysis of Tetramethylene Disulfotetramine in Foods Using Solid-Phase Microextraction–Gas Chromatography–Mass Spectrometry. *J. Chromatogr. A* **2008**, 1192, 36–40.
10. Saito-Shida, S.; Nemoto, S.; Teshima, R.; Akiyama, H. Determination of Rodenticide Tetramethylenedisulfotetramine (Tetramine) in Processed Foods by Gas Chromatography–Tandem Mass Spectrometry. *Shokuhin Eiseigaku Zasshi [Journal of the Food Hygienic Society of Japan]* **2016**, 57 (3), 72–75.
11. Owens, J.; Hok, S.; Alcaraz, A.; Koester, C. Quantitative Analysis of Tetramethylenedisulfotetramine (Tetramine) Spiked into Beverages by Liquid Chromatography–Tandem Mass Spectrometry with Validation by Gas Chromatography–Mass Spectrometry. *J. Agric. Food Chem.* **2009**, 57 (10), 4058–4067.
12. Arthur, C.L.; Pawliszyn, J. Solid Phase Microextraction with Thermal Desorption Using Fused Silica Optical Fibers. *Anal. Chem.* **1990**, 62 (19), 2145–2148.

Blank

## ACRONYMS AND ABBREVIATIONS

CRM	consecutive reaction monitoring
ESI	electrospray ionization
GC	gas chromatography
IT-MS	ion trap mass spectrometry
LC	liquid chromatography
LD <sub>50</sub>	lethal dose for 50% of test subjects
MS	mass spectrometry
<i>m/z</i>	mass-to-charge ratio
SRM	selected reaction monitoring
TEA	trimethylamine
TETS	tetramethylenedisulfotetramine
USDA	U.S. Department of Agriculture


Blank

## APPENDIX

### SELECTED REACTION MONITORING (SRM) AND CONSECUTIVE REACTION MONITORING (CRM) CHROMATOGRAMS AND MASS SPECTRA OF TETRAMETHYLENEDISULFOTETRAMINE (TETS) EXTRACTED FROM VARIOUS FOOD MATRICES

#### A.1. TETS and Apple Juice


RT: 0.00 - 4.99 SM: 15G

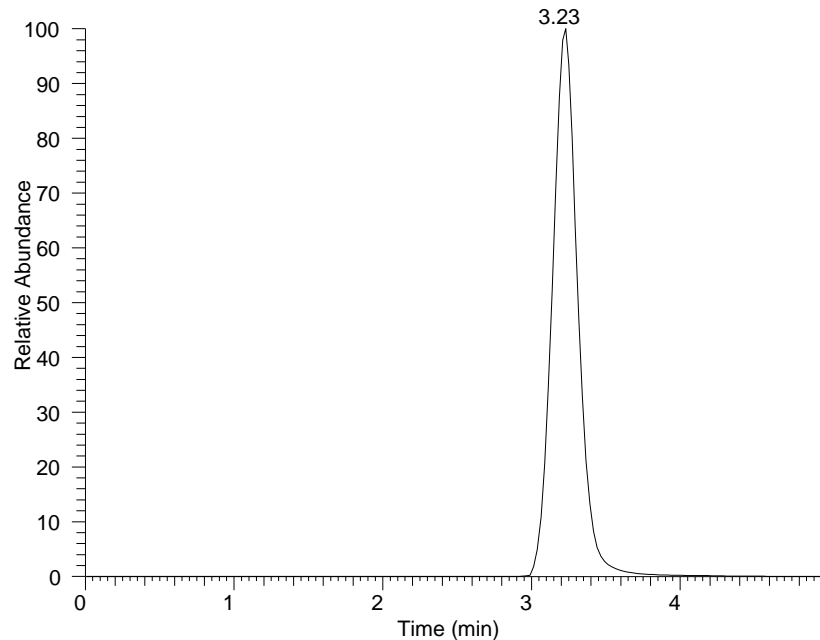

NL:  
7.90E4  
TIC F: ITMS - c  
ESI SRM ms2  
347.00@cid15.00  
[266.50-269.50]  
MS AJ1-1


AJ1-1 #250-272 RT: 3.11-3.34 AV: 11 NL: 5.77E4

F: ITMS - c ESI SRM ms2 347.00@cid15.00 [266.50-269.50]

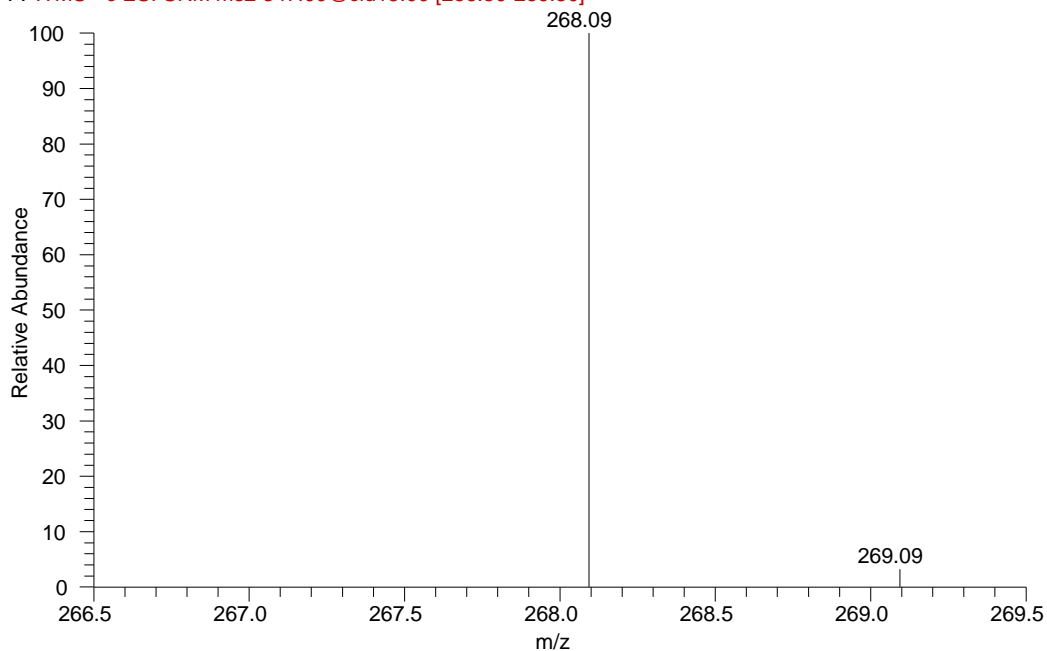

## A.2. TETS and Orange Juice


RT: 0.00 - 5.00 SM: 15G

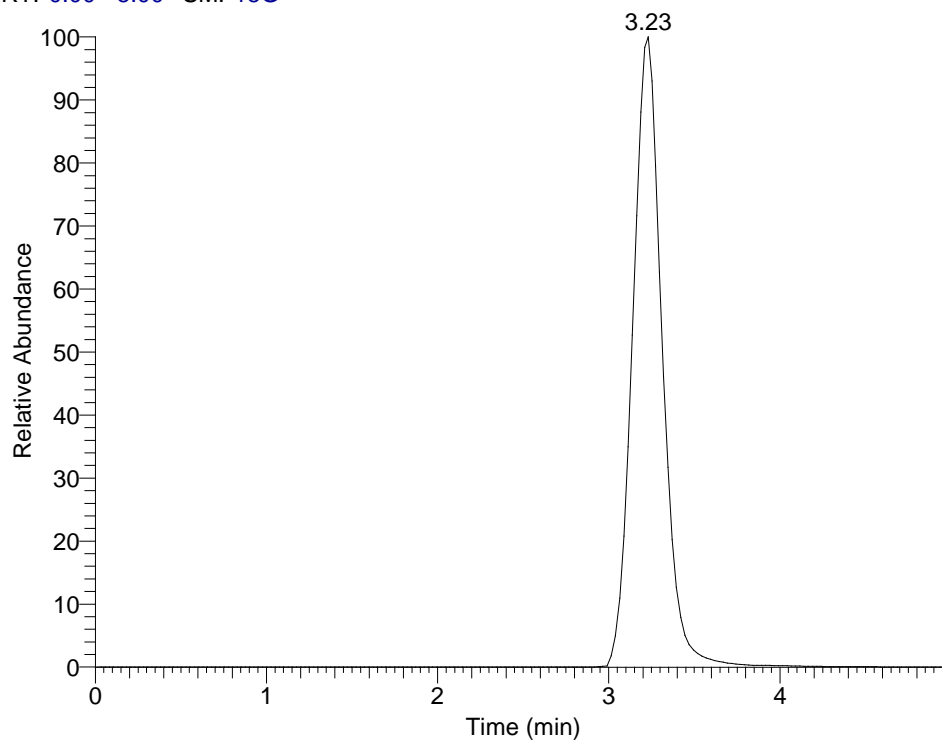

NL:  
6.73E4  
TIC F: ITMS - c  
ESI SRM ms2  
347.00@cid15.00  
[266.50-269.50]  
MS OJ1-1


OJ1-1 #251-271 RT: 3.11-3.35 AV: 11 NL: 4.90E4

F: ITMS - c ESI SRM ms2 347.00@cid15.00 [266.50-269.50]

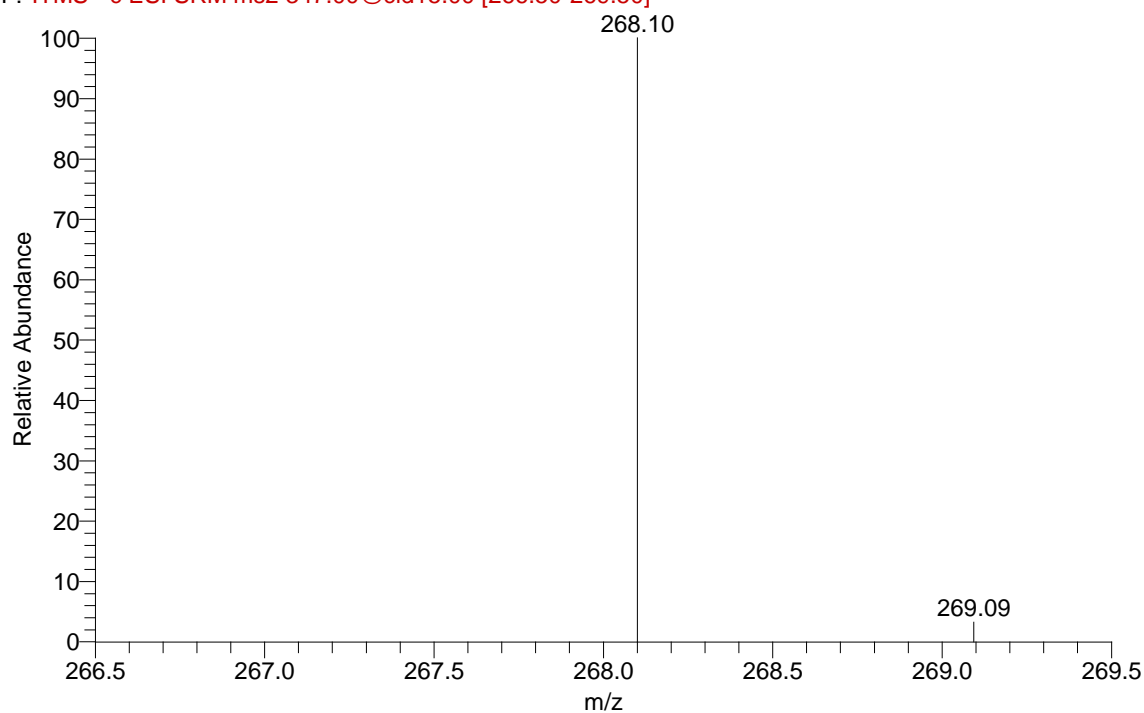


### A.3. TETS and Egg Beaters Egg Whites

RT: 0.00 - 5.00 SM: 15G

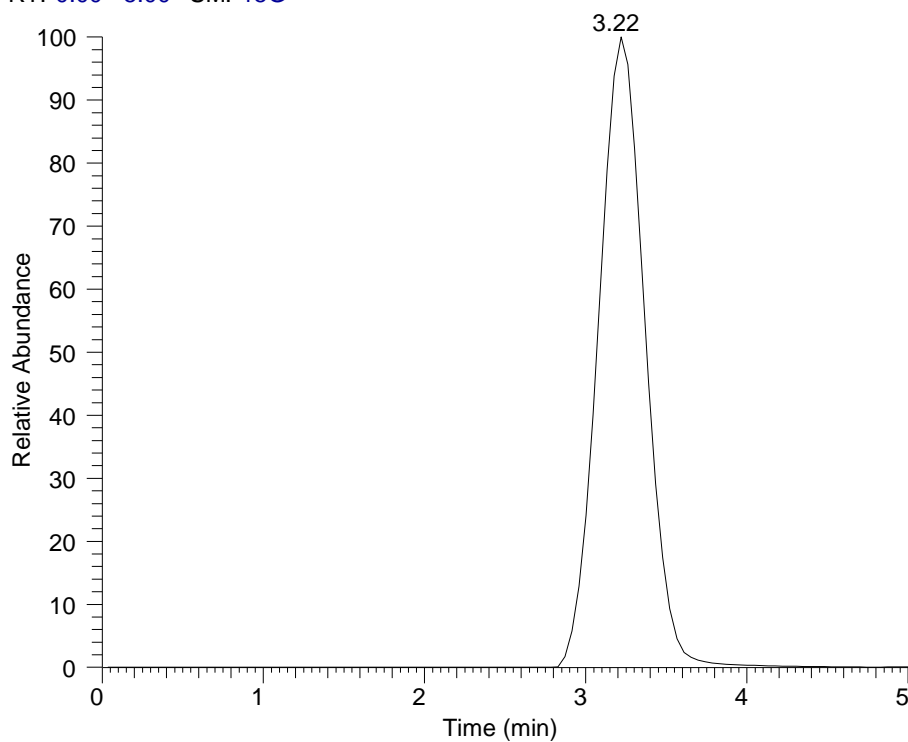


NL: 1.75E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
EB-1-1

EB-1-1 #213-233 RT: 3.09-3.35 AV: 7 NL: 1.48E4  
F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

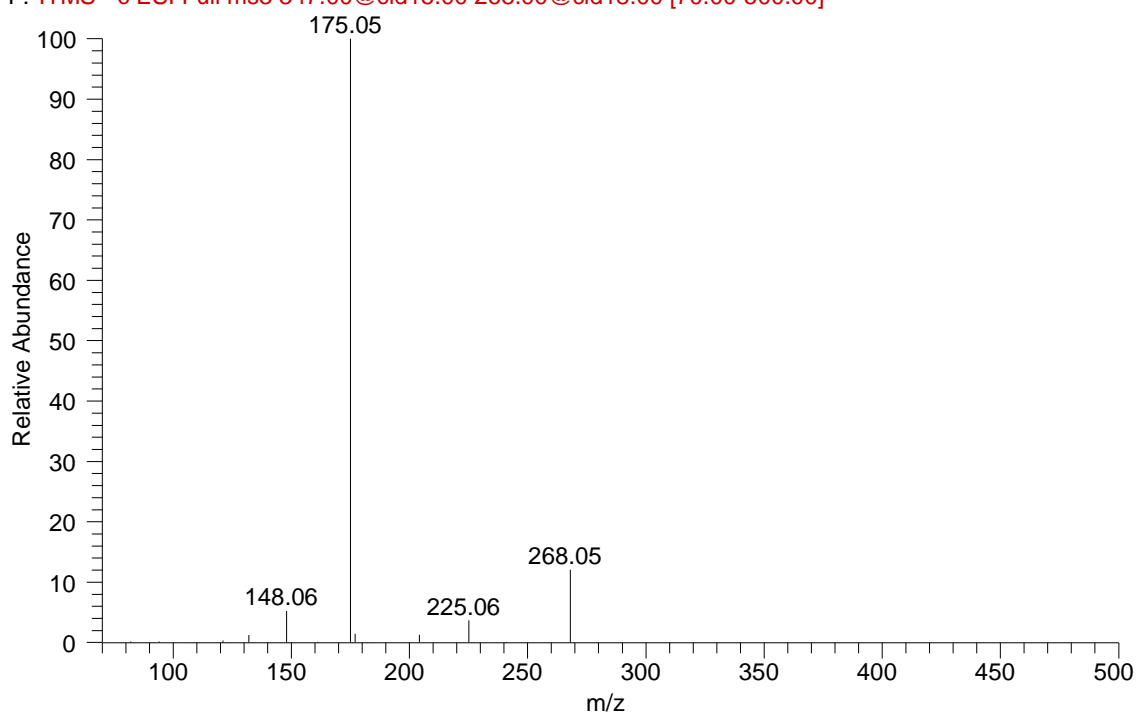


#### A.4. TETS and USDA Egg White

RT: 0.00 - 5.00 SM: 15G

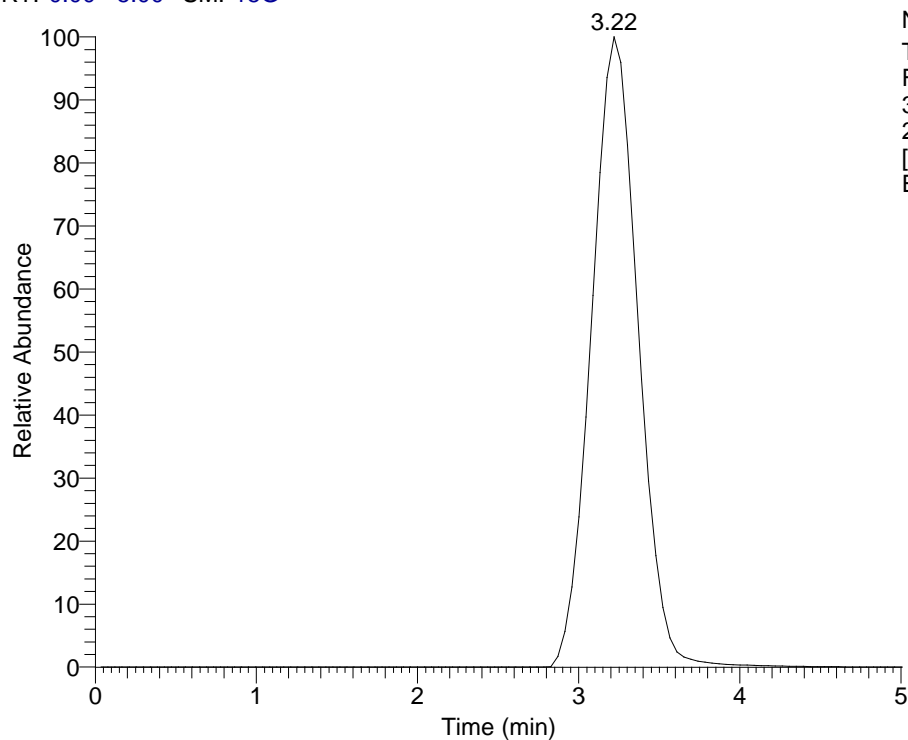


NL: 1.77E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
EW-1-1

EW-1-1 #209-236 RT: 3.05-3.39 AV: 9 NL: 1.18E4  
F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

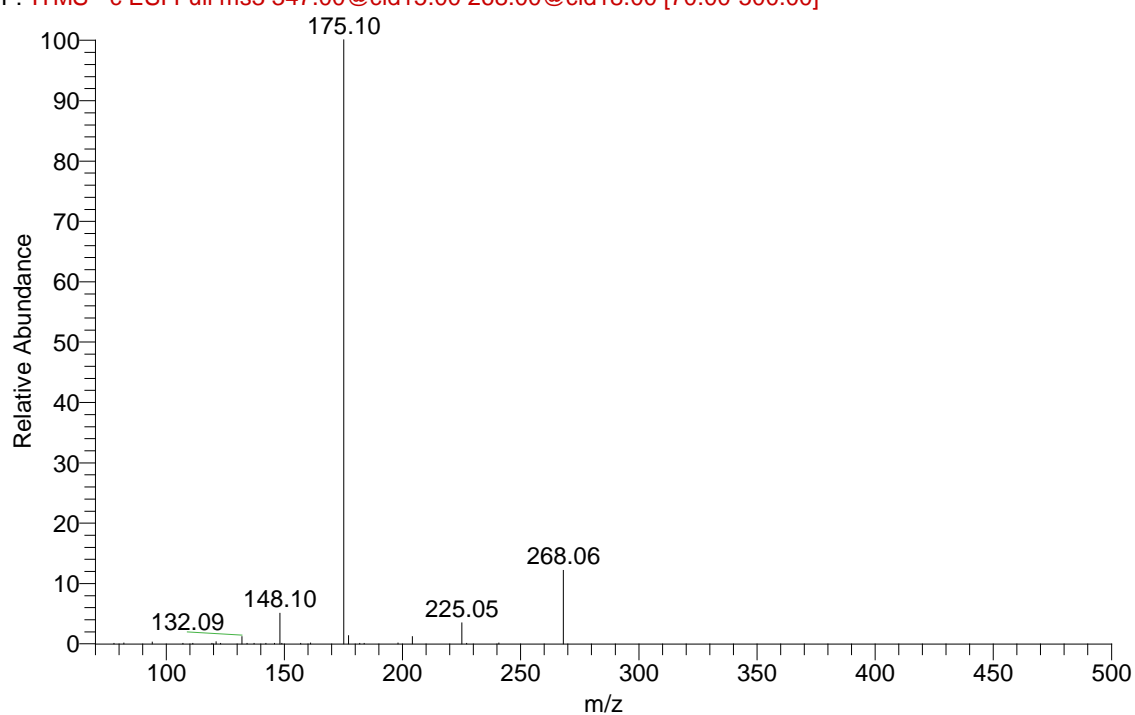

## A.5. TETS and USDA Egg Yolk 2551


RT: 0.00 - 5.00 SM: 15G

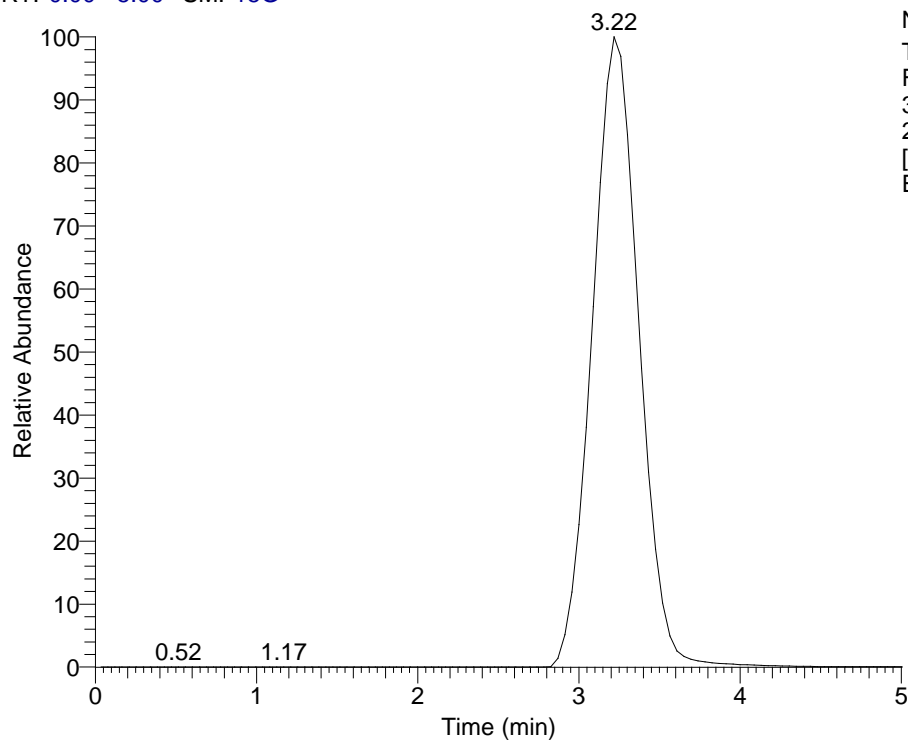

NL: 1.85E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
EY-1-1


EY-1-1 #212-235 RT: 3.09-3.39 AV: 8 NL: 1.39E4

F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

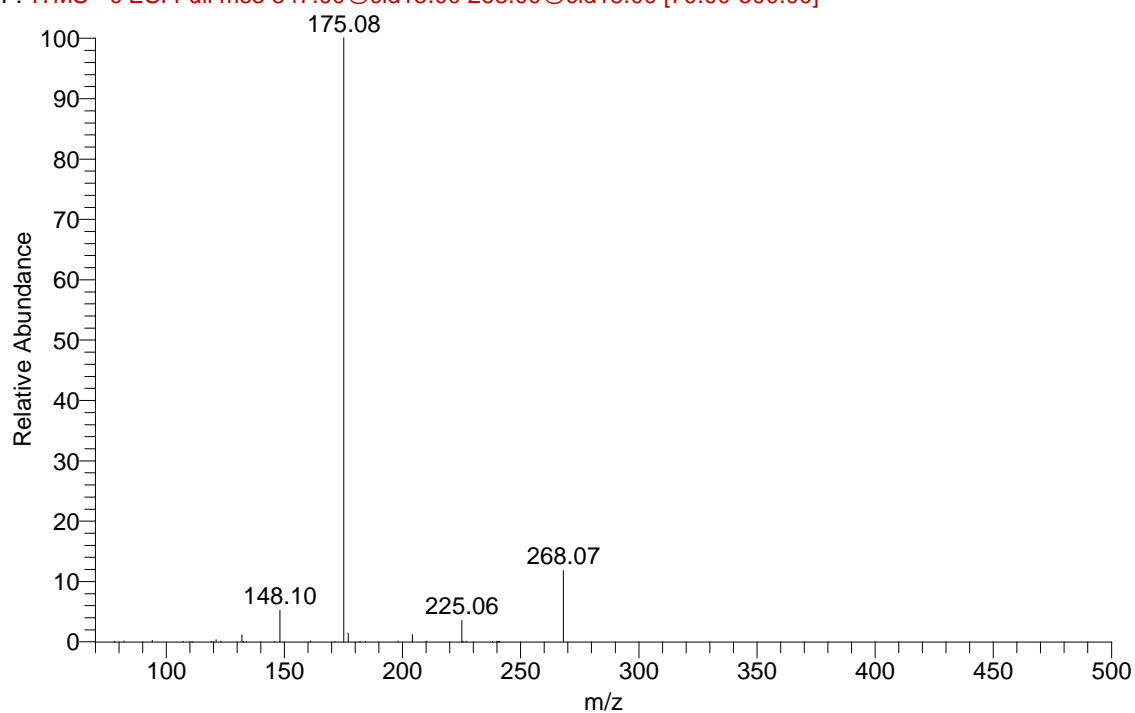

## A.6. TETS and USDA Egg Yolk C24410


RT: 0.00 - 5.00 SM: 15G

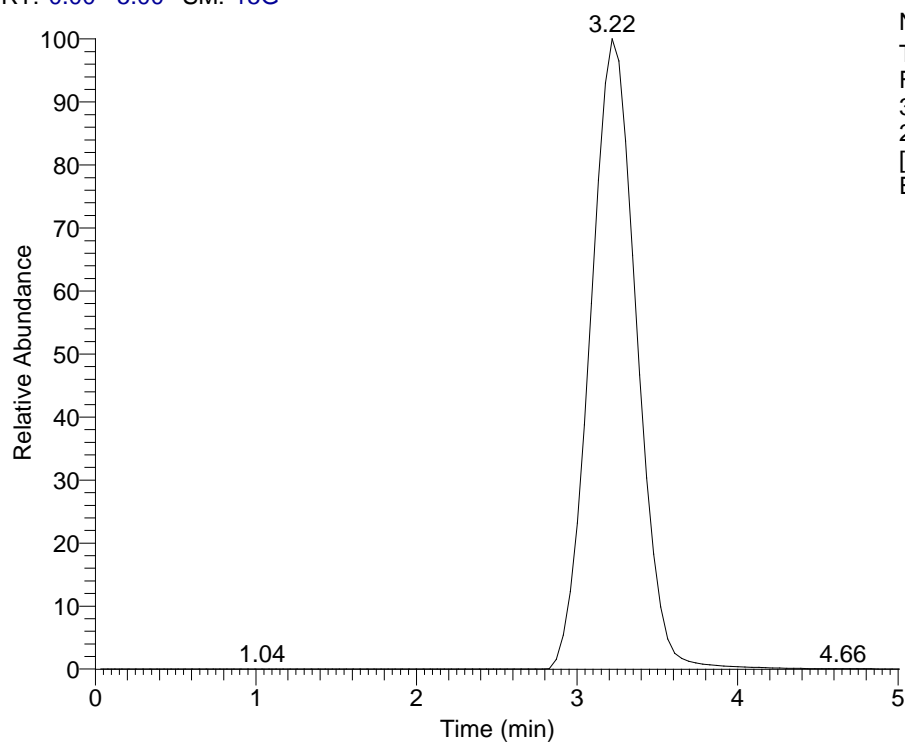

NL: 1.83E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
EYC-1-1


EYC-1-1 #214-231 RT: 3.13-3.35 AV: 6 NL: 1.81E4

F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

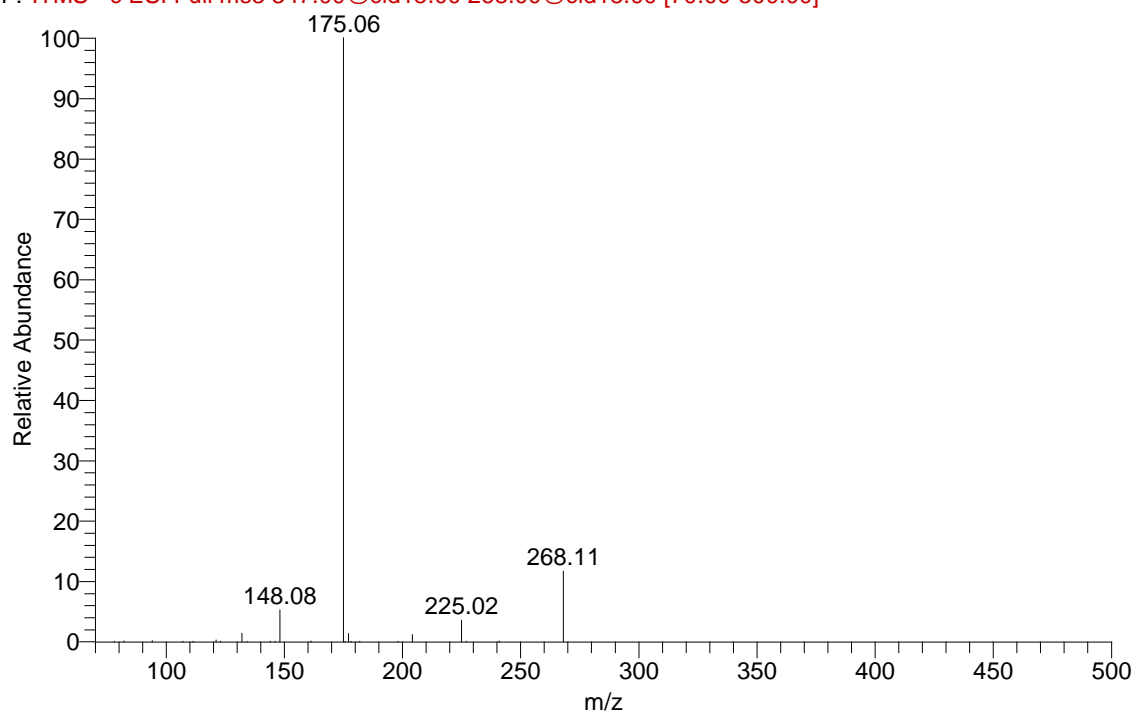

## A.7. TETS and Whole Egg


RT: 0.00 - 5.00 SM: 15G

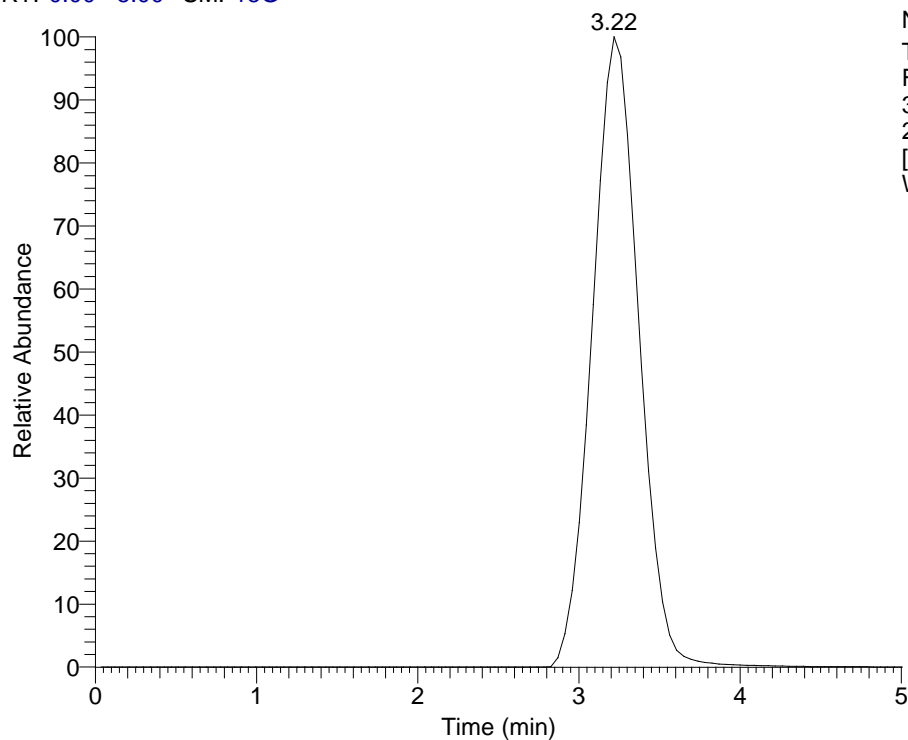

NL: 1.76E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
WE-1-1


WE-1-1 #210-236 RT: 3.05-3.39 AV: 9 NL: 1.19E4

F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

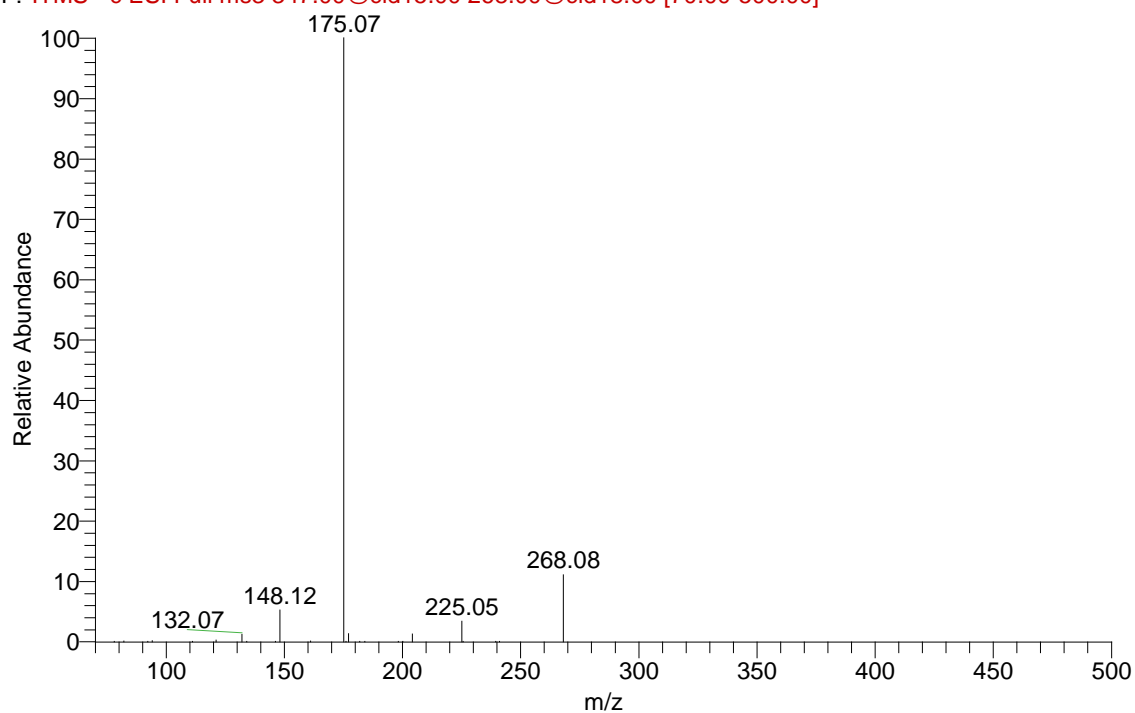

## A.8. TETS and Chicken Nuggets


RT: 0.00 - 5.00 SM: 15G

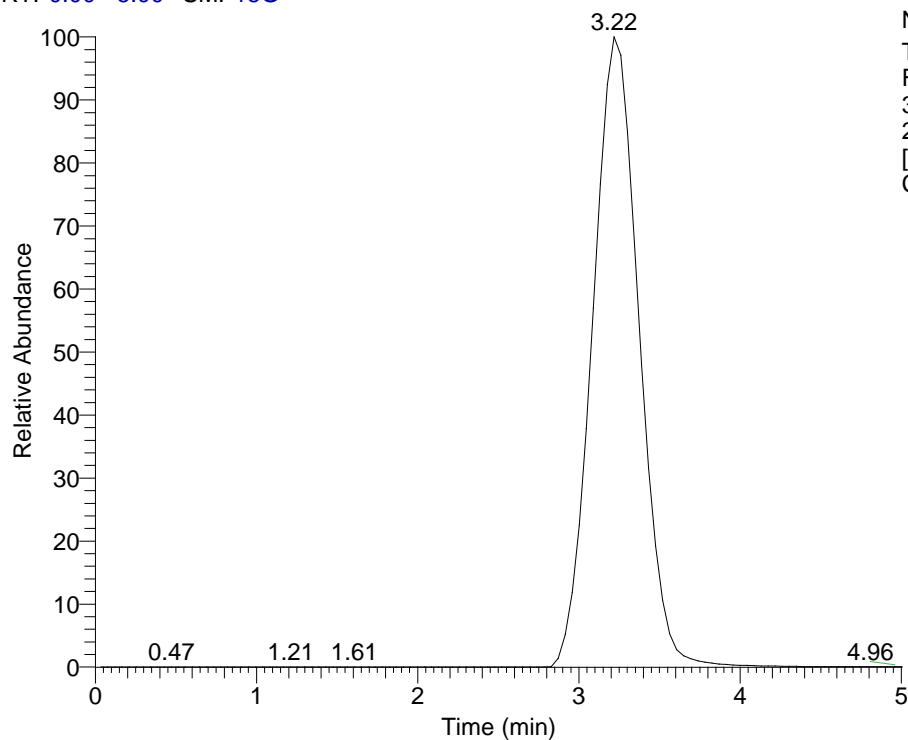

NL: 2.01E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
CN-1-1


CN-1-1 #207-240 RT: 3.00-3.48 AV: 12 NL: 1.03E4

F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

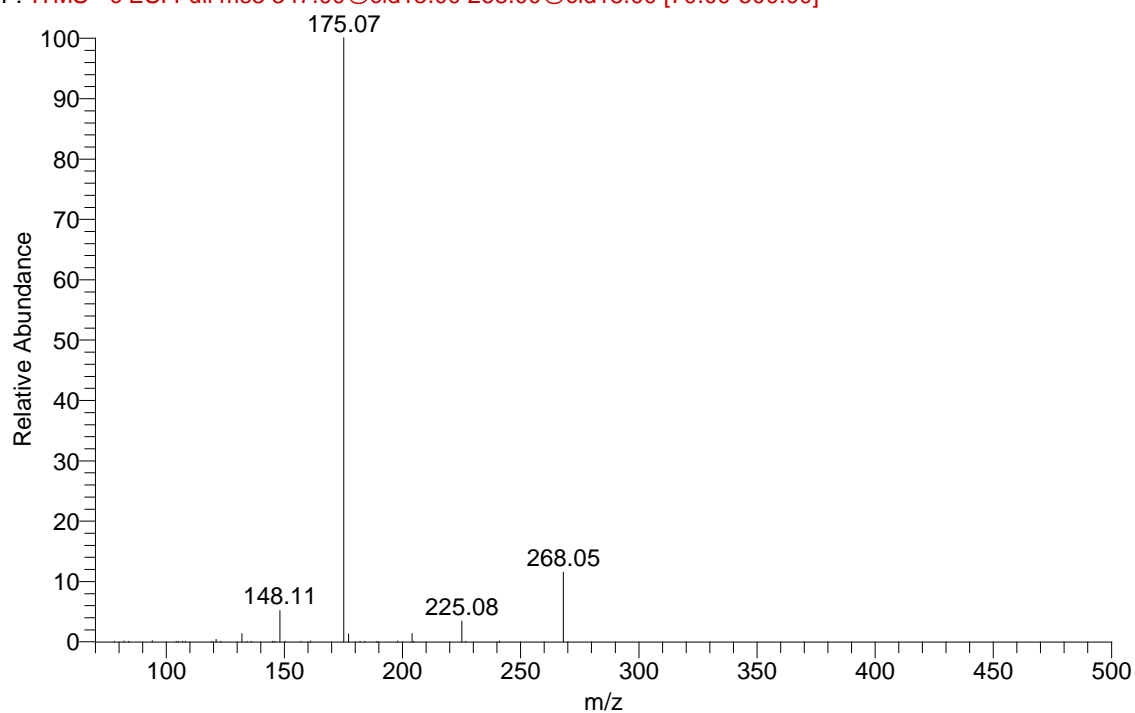


## A.9. TETS and 80/20 Ground Beef

RT: 0.00 - 5.00 SM: 15G

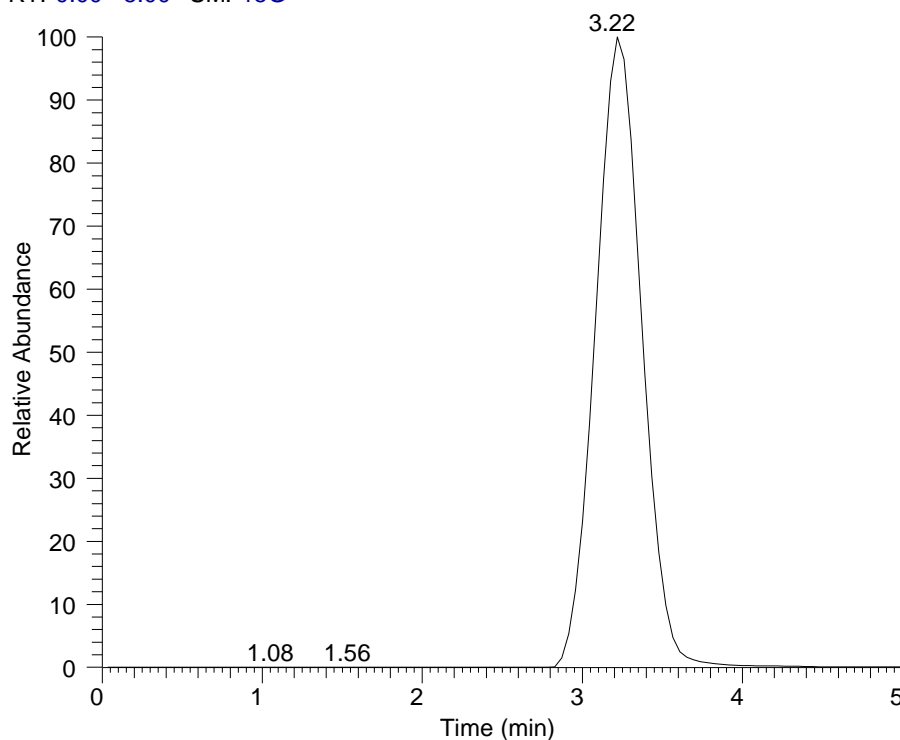


NL: 1.79E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
BEEF-1-1

BEEF-1-1 #212-235 RT: 3.09-3.39 AV: 8 NL: 1.35E4  
F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

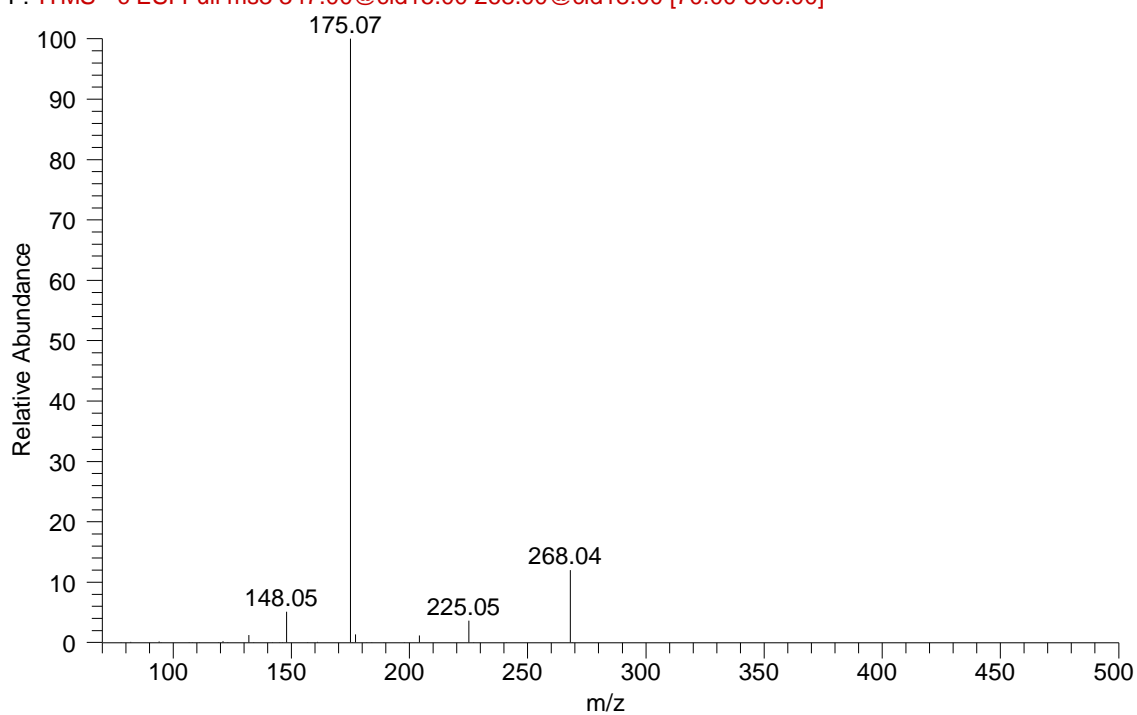

## A.10. TETS and Hot Dog


RT: 0.00 - 5.00 SM: 15G

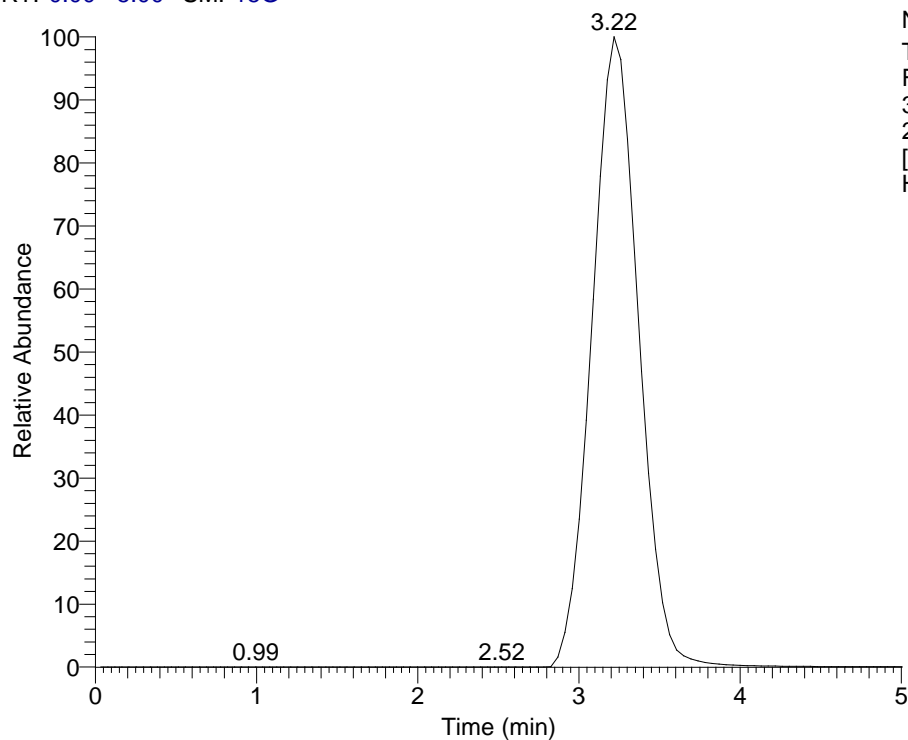

NL: 2.12E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
HD-1-1


HD-1-1 #212-235 RT: 3.09-3.39 AV: 8 NL: 1.60E4

F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

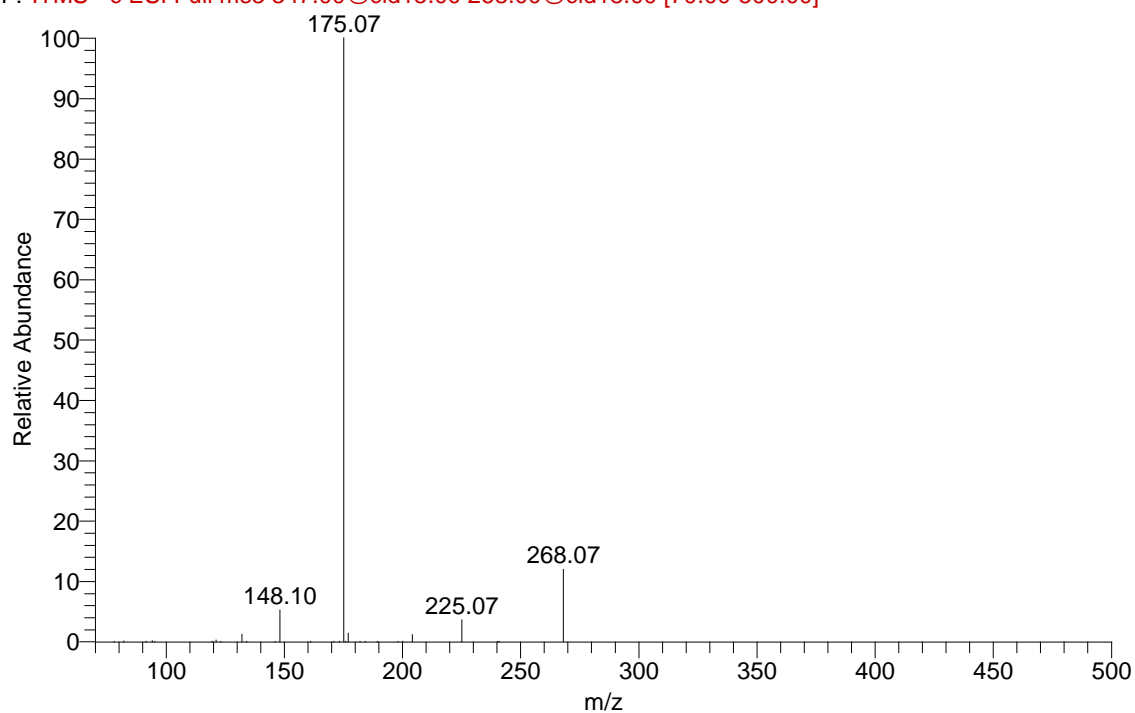

## A.11. TETS and Turkey Deli Meat


RT: 0.00 - 5.00 SM: 15G

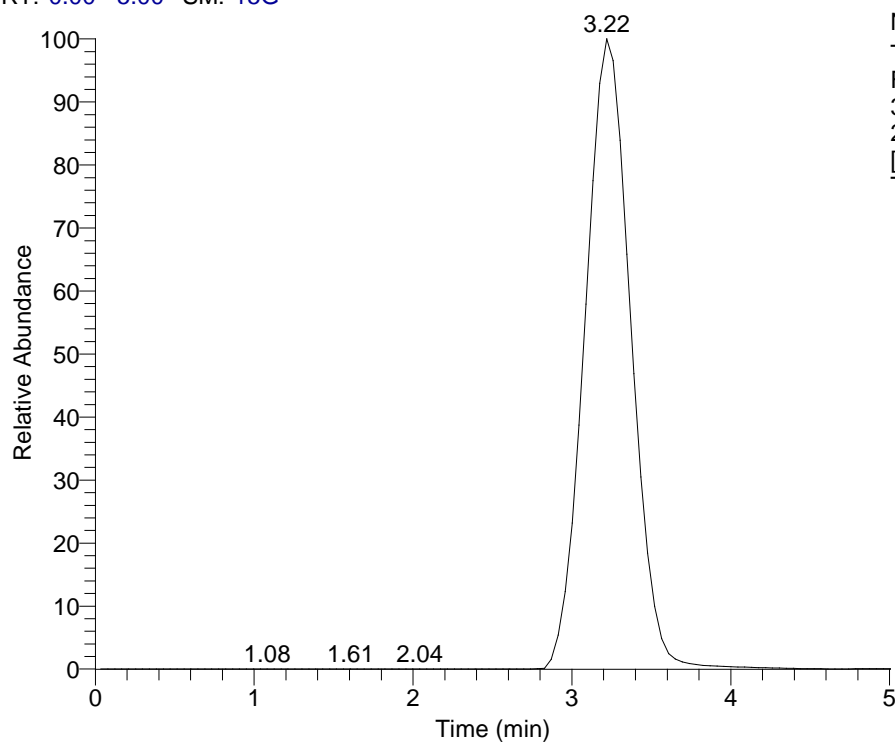

NL: 1.70E4  
TIC F: ITMS - c ESI  
Full ms3  
347.00@cid15.00  
268.00@cid18.00  
[70.00-500.00] MS  
TURKEY-1-1

TURKEY-1-1 #212-235 RT: 3.09-3.39 AV: 8 NL: 1.29E4

F: ITMS - c ESI Full ms3 347.00@cid15.00 268.00@cid18.00 [70.00-500.00]

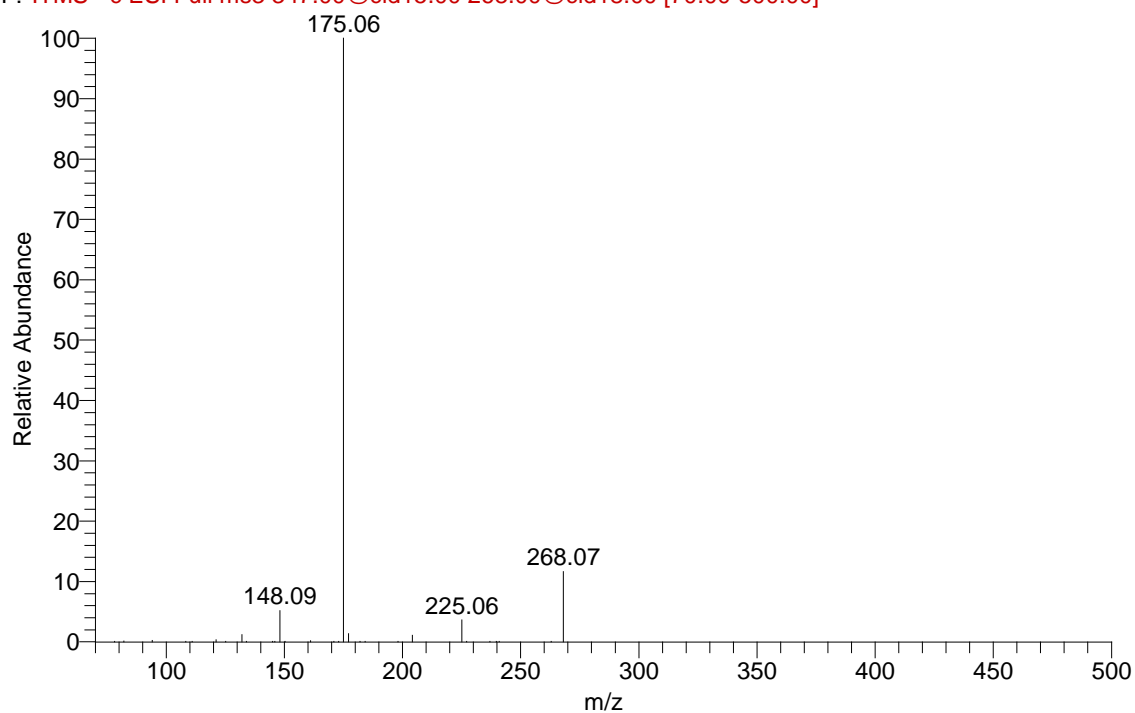



## **DISTRIBUTION LIST**

The following individuals and organizations were provided with one Adobe portable document format (pdf) electronic version of this report:

U.S. Army Edgewood Chemical  
Biological Center (ECBC)  
RDCB-DRC-C  
ATTN: Bae, S.  
Kong, L.  
McMahon, L.

ECBC Technical Library  
RDCB-DRB-BL  
ATTN: Foppiano, S.  
Stein, J.

Defense Technical Information Center  
ATTN: DTIC OA

Defense Threat Reduction Agency  
J9-CBS  
ATTN: Graziano, A.

Department of Homeland Security  
RDCB-PI-CSAC  
ATTN: Negron, A.  
DHS-S&T-RDP-CSAC  
ATTN: Strang, P.

G-3 History Office  
U.S. Army RDECOM  
ATTN: Smart, J.

Office of the Chief Counsel  
AMSRD-CC  
ATTN: Upchurch, V.

ECBC Rock Island  
RDCB-DES  
ATTN: Lee, K.  
RDCB-DEM  
ATTN: Grodecki, J.

